

Functional Elbow Brace

Hand-built by Thuasne USA, the Functional Elbow brace is raising the industry standard for protection and suspension with a brace engineered to provide ultimate protection for the unique anatomy of the elbow.





Functional Elbow


Recommended for functional prophylactic use. Chronic elbow injuries, mild to moderate soft tissue injuries, the treatment of stable fractures, or for postoperative use.

Thuasne USA's Functional Elbow

combines trusted hinge technology with a three-dimensional uni-frame structure for strength along with moldable metal cuffs, and Thuasne's patented Compression/ Suspension (C/S) technology.

- Universal left/right
- Anti-migration/anti-microbial pads
- Adjustable range-of-motion (ROM)
 Reversible strapping pattern
- Removable medial condylar suspension wedge

C/STECHNOLOGY

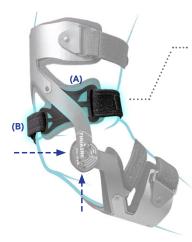

Compression/ Suspension (C/S)

Anti-migration technology made up of a semi-rigid padded Sub-shell (A) and adjustable strap (B), designed to maintain proper hinge alignment.

5545 St. Jacques O. Montreal, QC H4A 2E3 O: 514. 369. 3311 Toll Free / Sans Frais : 1.877.527.1911 /Fax Toll Free Telec Sans Frais : 1.877.527.3311 www.bmiortho.com

Strength • Suspension • Fit

Working in unison with athletic trainers, clinicians and Thuasne's team of experts, the Functional Elbow is built from the ground up to halt migration, improve suspension and protect while providing a customized fit perfect for athletes and patients.

FEATURES



*Elbow circumference is best taken in full extension at the medial epicondyle



SIZE	ELBOW CIRCUMFERENCE
S	9.75 - 11.25
М	11.25 - 12.75
L	12.75 - 14.25
XL	14.25 - (+)

If a patient has more soft tissue, recommendation is sizing down if they are within .25-.5" of the lower end of the circumference measurement.