



**THUASNE**

# SafeLimb

Trans tibial Residual Limb Protector

**BMI**  
BERRETTA MEDICAL INC.



## Features

Post-surgical  
support & protection

Five sizes; universal  
for left or right leg

Bars can be  
hand-contoured &  
lengthened

Anterior opening  
& straps for easy  
application

Cushioned rigid plate  
under residual limb

# SafeLimb


## Transtibial Residual Limb Protector

### RECOMMENDED L5450


During post-surgical care of transtibial amputation patients, the SafeLimb provides supportive protection of the residual limb if the patient falls.

Prosthetists will appreciate the SafeLimb's unique design, materials and adaptability - which enhances the fit, suspension and patient compliance. A spherical compressive foam cushion mounted on a rigid plate - attached to the struts at the distal end of the brace - prevents trauma to the surgical wound.

The length of the SafeLimb's medial and lateral bars is adjustable, and the bars can also be hand-contoured to align with the degree of extension. The material interface is lightweight and breathable, with anterior-opening and anterior straps for ease of application. A proprietary "SupraPatellar Strap" (above the knee) effectively prevents migration.

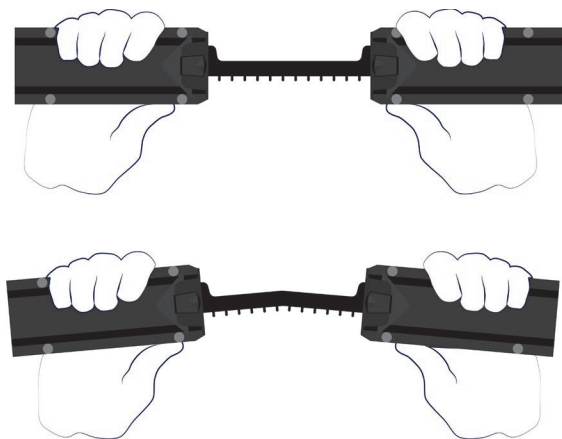

*The SafeLimb's metal bars can be contoured without tools.*

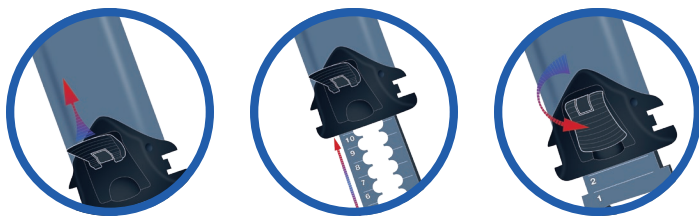

*Townsend's proprietary "SnapLock" allows for quick, tool-free length adjustment.*



*The residual limb is protected by a cushioned plate at the distal end of the brace*

Size	Thigh Circumference 6" Above Mid Patella
SM	14 – 18 ¾"
MD	18 ¾ – 20 ¾"
LG	20 ¾ – 24 ¾"
XL	24 ¾ – 29 ½"
XXL	29 – 34"



5545 St. Jacques O. Montreal, QC H4A 2E3  
O: 514. 369. 3311

Toll Free / Sans Frais : 1.877.222.3311 ✶ Fax Toll Free / Telec Sans  
Frais : 1.877.527.1911 [www.bmiortho.com](http://www.bmiortho.com)